



***KiaTo*** 

The smart choice for sharp blades



Kiato Plus scalpels are made with passion using a state-of-the-art robotized production line in compliance with all regulatory requirements.

Kiato Plus scalpels set the benchmark for all others to follow, from East to West. Rich on features and high on performance, Kiato Plus scalpels are user-friendly with a laser controlled durable sharp edge. In spite of its premium features, use of these scalpels enables hospitals to scale down their budgets. It is, therefore, not surprising that the Kiato Plus range of products is establishing a strong presence in markets worldwide.

Kiato Plus disposable scalpels are manufactured in an ISO 13485:2016 certified factory and are CE marked and gamma sterilised. They are individually packed, 10 to a box, in easy peel pouches. Each of these scalpels have a premium Swedish stainless steel surgical blade fitted to a plastic handle with a safety guard for safe handling.

User-friendly  
Ergonomic and safe  
Premium quality  
Superior performance



REF. 23110




**KIATO** 



## User-friendly

The range of Kiato Plus scalpels come in easily differentiated color-coded boxes. Each of these fibre-free, easy peel packs display the blade shape for ease of identification. The blade size is also prominently printed on the handle with a graduated cm scale. Ease of use begins with the product packaging itself and makes for greater convenience.

REF. 23111

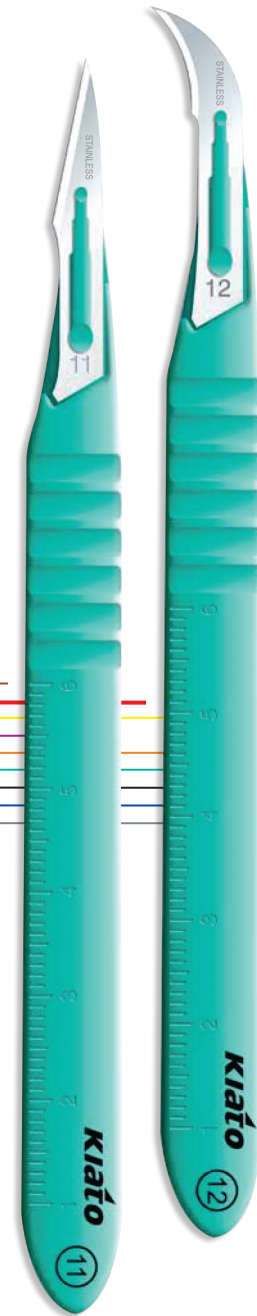

REF. 23112



## Ergonomic and safe

Made of ABS, the Kiato Plus scalpel handles are not just ergonomically designed but also stiffer and heavier, to facilitate ease of use. A serrated grip ensures greater comfort and control while the scale on the handle helps in determining incision length. For enhanced safety, each scalpel also has a unique stiff cap with lock tabs that prevents injuries and improves performance.



REF. 23115




REF. 23120

15

20

REF. 23121

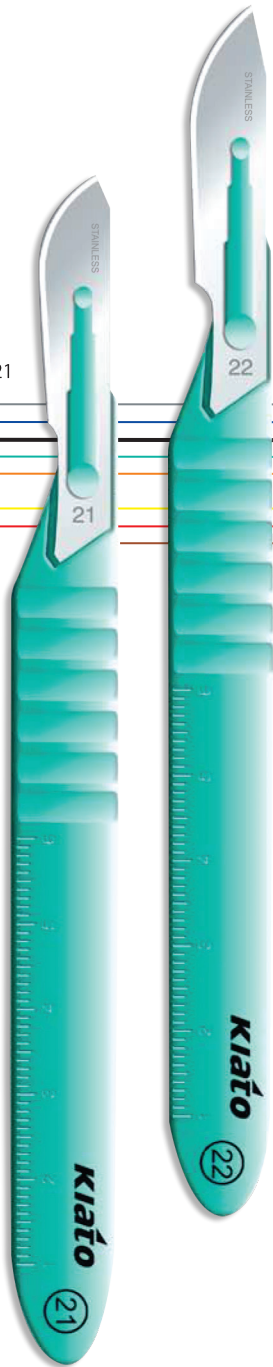

REF. 23122

## Premium quality

The scalpel handles are molded in a full hot runner mold. The robotized handling system prevents any damage to the handle by gently placing it on the conveyor. Our agreement with Sandvik, Sweden ensures that we receive and use superior stainless steel quality. The blades are made of Sandvik grade 13C26 which conforms to the undermentioned specifications.

Stainless Steel Chemical Composition  
Ref: BS 5194: 1991, Type F

Constituents	Percentage (%)
Carbon (C)	0.60 to 0.70
Silicon (Si)	0.50 maximum
Manganese (Mn)	1.00 maximum
Chromium (Cr)	12.00 to 13.50
Phosphorous (P)	0.03 maximum
Sulphur (S)	0.025 maximum
Nickel (Ni)	0.50 maximum




**kiato** 



REF. 23123  
10 Sterile Einmalskalpelle  
10 Bistouris stériles à usage unique  
10 Bisturis estériles para uso unico  
10 Bisturi sterili monouso  
10 Sterila engångsskalpeller

Stainless Steel Surgical blade fitted on a plastic handle with safety shield

REF. 23123

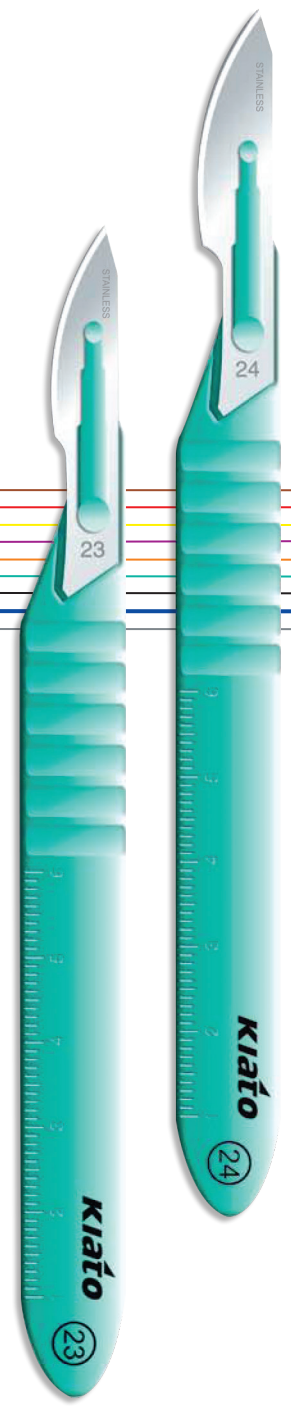

REF. 23124



## Superior performance

Sharpness and durability are the true tests of a good quality scalpel. Both these factors are hallmarks of every Kiato Plus scalpel. The laser controlled edge in every scalpel guarantees benchmark sharpness and durability to ensure superior customer satisfaction and unflinching performance, even over time.



Aditya Dispomed Products Pvt Ltd  
Plot No. 19 Sector 6 I.M.T. Manesar  
Gurgaon Haryana 122 050 India  
[www.adityadispomed.com](http://www.adityadispomed.com)



MDSS GmbH  
Schiffgraben 41  
30175 Hannover  
Germany

